




THAI FURNACES ENGINEERING

HI-TEMP

## Digital Convection Dry Oven

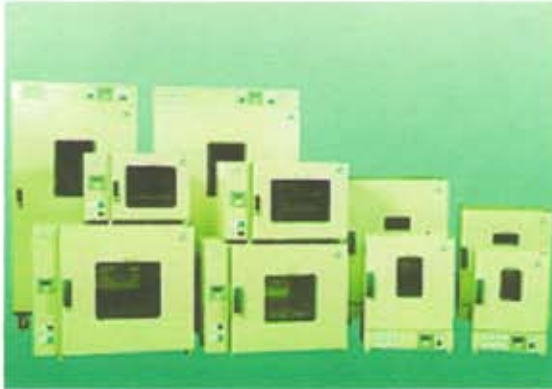

### 1 Features

- Designed for general laboratory use, these ovens feature a double wall for durability and excellent heat retention.
- The high temp. silicone rubber gasket prevent heat loss.
- Temperature is controlled by Digital PID
- The oven features a dial with increments for set point temperature.
- The stainless steel shelves are adjustable, and view glass let you the samples during operation.



### 2 Specifications

Model		HV-70L	HV-135L	HV-220L
Capacity		72 Liter	135 Liter	220 Liter
Chamber WxLxH(mm.)		400x400x450	450x550x550	500x600x750
Temp.	Range	Ambient temperature +10 to 300°C		
	Control	Digital PID control with Pt 100 ohm sensor		
Temperature fluctuation		less than 0.5°C or better		
Heater		<1400W	<2000W	<2600W
Material	Interior	Stainless steel		
	Exterior	Lacquer coated on high quality steel plate		
	Insulation	Superior glass fiber		
Timer		9999 minutes		
Air circulation		with a fan motor		
Control & Display		Microprocessor digital PID control & LCD display		
Door		Silicon gasket between the chamber and door		
Electric source		220VAC, 50/60Hz		

สอบถามรายละเอียดเพิ่มเติมได้ที่ 087-7117701 คุณมนตรี, 087-1770404 คุณทศพร